

Cuvette Test **LCK 381**

TOC Total organic carbon

Total carbon (TC) and total inorganic carbon (TIC) are converted to carbon dioxide (CO₂) by, respectively, oxidation and acidification. The CO₂ passes from the digestion cuvette through a membrane and into the indicator cuvette. The change of colour of the indicator is photometrically evaluated. TOC (total organic carbon) is determined as the difference between the TC and TIC values.

Range of Application

Waste water, surface water, soils

Storage Information

The test reagents are stable at +15 to +25°C up to the expiry date given on the package.

Special note

- 1. Contamination by ambient air
 - Never leave cuvette open, because carbon dioxide in the ambient air can cause high bias. Cuvettes must only be opened when necessary (e.g. to add sample) and must be closed or further processed immediately afterwards.
- 2. How to use the powder dispenser Screw powder dispenser on to digestion reagent A (LCK 381 A). Invert so that powder dispenser is under the reagent and shake. This causes the dispensing chamber to be filled. Position the centring recess of the powder dispenser above the TC cuvette and add one dose. Close digestion reagent A (LCK 381 A) immediately with original cap.
- 3. Labelling the digestion cuvettes If several samples are analysed simultaneously, label them so that the TC and TIC cuvette combinations of the same sample
- 4. Thermostat

Heat the thermostat to 100°C (check the temperature - higher temperatures lead to dangerous excess pressure). When this temperature has been reached, insert the cuvette combinations and start the reaction time (2h) again. Insert cuvette combinations only in the small shafts in the thermostat. Do not insert them in the large shaft with reducing sleeves.

5. Disposal

Do not screw the cuvette combinations apart when the analysis has been completed, but press them back into the blister pack in combination with the indicator cuvette upwards.

Interferences

If the TC and/or the TIC results are higher than the upper limit of the measuring range, the calculated TOC result which is displayed may be within the measuring range.

In this case and generally the measurement results must be subjected to plausibility checks (dilute and/or spike the sample).

Use only carbon-free double-distilled water to dilute the sample. The ions listed in the table have been individually checked up to the given concentrations. Cumulative effects and the influence of other ions have not been determined by us.

There is no interference from:

TIC-determination:

800 mg/L: HCOO 500 mg/L: CH₂COO 60 mg/L: SO₃2-20 mg/L: S2-6 mg/L: NO₂-N

Higher concentrations of these ions cause high-bias results.

TC-determination:

5000 mg/L: Cl 2000 mg/L: Ca²⁺, Mg²⁺ 1000 mg/L: NH₄-N

Higher concentrations of these ions cause low-bias results.

pH/Temperature

The pH of the sample must be between pH 4 and pH 10. The temperature of the sample and reagents must be between 15 and 25°C.

Safety Advice

On grounds of quality and reliability, the analysis should be carried out only with original HACH LANGE accessories.

Be sure to set the required temperature to 100°C (at 148°C the cuvette combinations may break apart).

The digestion conditions are such that oxygen is formed, and this causes a rise in pressure in the cuvette combinations. If the cuvette combinations are subjected to strong mechanical stress after the digestion reaction, e.g. if they suffer a blow or a fall, they may shatter. In this case glass splinters may cause injury.

CADAS 100 (LPG 185 / ≥ LPG 210)

If this test is not already stored in your instrument please ask your HACH LANGE Agency for programming instructions.

Procedure LCK 381

Applies to all types of photometer

TOC Total organic carbon

Edition 96/03

Please read carefully the information under the heading "Special note".

Carry out each TC determination in sequence in order to avoid contamination by ambient air (diagram A – E).

1. TC determination (see diagram A - E)

Into the TC cuvette dispense and pipette:

Digestion reagent A (LCK 381 A)

1 dose of powder

Sample

0.2 n

Close **TC cuvette** with original cap and invert a few times. Close the indicator cuvette **very tightly** with the membrane double-cap (screw on the membrane double-cap in such a way that the barcode label is on the lower half of the cap). **Immediately** screw the prepared indicator cuvette **tightly** on to the **TC cuvette**.

NB: Cuvettes must be held vertically and must not be inverted! The sample must not come into contact with the membrane of the double-cap!

Procedure LCK 381

Applies to all types of photometer

TOC Total organic carbon

Edition 96/03

1 mL

Please read carefully the information under the heading "Special note".

Carry out each TIC determination in sequence in order to avoid contamination by ambient air (diagram F – I).

2. TIC determination (see diagram F - I)

Into the TIC cuvette pipette:

Sample

Close **TIC** cuvette with original cap and invert a few times. Close the indicator cuvette *very tightly* with the membrane double-cap (screw on

indicator cuvette **very tightly** with the membrane double-cap (screw on the membrane double-cap in such a way that the barcode label is on the lower half of the cap). **Immediately** screw the prepared indicator cuvette **tightly** on to the **TIC cuvette**.

NB: Cuvettes must be held vertically and must not be inverted! The sample must not come into contact with the membrane of the double-cap!

3. Heat (see diagram J + K)

Place the cuvette combinations in the *preheated* thermostat with the coloured indicator cuvette upwards, and heat for *2 h* at *100°C*. For **TOC** determination the **TC** and **TIC** cuvette combinations must be heated simultaneously in the thermostat. Allow to cool to room temperature (indicator cuvette upwards).

4. Evaluation (see diagram L + M)

First clean the outside of the indicator cuvette of **TC** cuvette combination thoroughly (NB: when you clean the cuvette, be careful not to turn it in such a way that the cuvette combination could accidentally be screwed apart). Place in photometer and measure – 1. Measurement. Then clean the outside of the indicator cuvette of **TIC** cuvette combination thoroughly (NB: when you clean the cuvette, be careful not to turn it in such a way that the cuvette combination could accidentally be screwed apart). Place in photometer and measure – 2. Measurement. Press the used cuvette combinations back into the blister pack (do not unscrew cuvette combinations!).

Data table LCK 381

LP2W 96/03

TC • $F_1 = 0$ • $F_2 = 1197$ • K = -245.6**TIC** • $F_1 = 0$ • $F_2 = 262.8$ • K = -45.94

CADAS 30/50 96/03

TC/TIC • λ: 435 nm • Pro.: 11 • F₁ = 1176 • F₂ = 0 • F₃ = 256.9 • F₄ = 0 •

 $K_1 = -283.1 \bullet K_2 = -53.05$

CADAS 30S/50S 96/03

TC/TIC • λ : 435 nm • Pro.: 11 • $F_1 = 1178$ • $F_2 = 0$ • $F_3 = 256.9$ • $F_4 = 0$ •

 $K_1 = -283.2 \bullet K_2 = -53.05$

ISIS 6000/9000 96/03

TC/TIC • λ : 430 nm • Pro.: 11 • F₁ = 1193 • F₂ = 0 • F₃ = 259.4 • F₄ = 0 •

 $K_1 = -298.1 \bullet K_2 = -56.43$

CADAS 100 / LPG 185 96/03

TOC/TC/TIC • λ : 435 nm • F_1 = 1165 • F_2 = -238.9 • F_3 = 253.8 •

 $F_4 = -43.16 \bullet F_5 = 1 \bullet F_6 = -1$

CADAS 100 / ≥ LPG 210 96/03

TOC/TC/TIC • λ : 435 nm • $F_1 = 1165$ • $F_2 = -238.9$ • $F_3 = 253.8$ • $F_4 = -43.16$

LASA 1 / plus

LCK 381

TOC Total organic carbon

Edition 96/03

Evaluation

- Press "Mode" key and check program control number:
 : 34
- 2. Insert program filter 440 nm.
- 3. Select test with "Mode" key.
- 4. Insert indicator cuvette of the **TC** cuvette combination.
- 5. Insert indicator cuvette of the **TIC** cuvette combination.

20 sec display = TC

20 sec display = TIC

20 sec display = **TOC**

or by pressing the key * the results for TC, TIC and TOC can be called up immediately one after another.

Parameter	Display	Meas. range
TOC	TOC LCK 381	60 – 735 mg/L
TC		75 – 750 mg/L
TIC		15 – 150 mg/L

LASA 10 / 20 **LCK 381**

TOC Total organic carbon

Edition 96/03

Evaluation

- 1. Press any key.
- 2. Check program control number: __: 34
- 3. Select test with \uparrow or \downarrow key.
- 4. Insert indicator cuvette of the **TC** cuvette combination.
- 5. Insert indicator cuvette of the **TIC** cuvette combination.

20 sec display = TC

20 sec display = TIC

20 sec display = TOC

or by pressing the key \downarrow the results for **TC**, **TIC** and **TOC** can be called up immediately one after another.

Parameter	Display	Meas. range
TOC	TOC LCK 381	60 – 735 mg/L
TC		75 – 750 mg/L
TIC		15 – 150 mg/L

LASA 30 **LCK 381**

TOC Total organic carbon

Edition 96/03

Evaluation

- 1. Insert filter 440 nm.
- 2. Select »Dr. Lange« mode.
- 3. Select test number (see below).
- 4. Control number must be 2.
- 5. Insert indicator cuvette of the TC cuvette combination and press green key.
- 6. Insert indicator cuvette of the TIC cuvette combination and press green key.

The result is displayed in TC, TIC and TOC.

Parameter	Test-No.	Meas. range
TOC	381	60 – 735 mg/L
TC		75 – 750 mg/L
TIC		15 – 150 mg/L

LP2W **LCK 381**

TOC Total organic carbon

Edition 96/03

Evaluation

TC determination

- 1. Insert program filter 435 nm.
- 2. Press "Tests" key until display (see below) appears.
- 3. Control number must be 8.
- 4. Insert blank-value cuvette (distilled water) and press "Null" (zero) key.
- 5. Insert indicator cuvette of the TC cuvette combination and press "Ergebnis" (result) key. Make a note of the display - TC

TIC determination


- 1. Insert program filter 435 nm.
- 2. Press "Tests" key until display (see below) appears.
- 3. Control number must be 4.
- 4. Insert blank-value cuvette (distilled water) and press "Null" (zero) key.
- 5. Insert indicator cuvette of the TIC cuvette combination and press "Ergebnis" (result) key. Make a note of the display - TIC

Calculation of the TOC concentration:

TC - TIC = mg/L TOC

Parameter	Display	Meas. range
TOC		60 – 735 mg/L
TC	TC LCK 381	75 – 750 mg/L
TIC	TIC LCK 381	15 – 150 mg/L



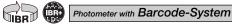

LCK 381

TOC Total organic carbon

Edition 96/03

Evaluation

- 1. Check program control number:
 - __: 34 (CADAS 30S/50S, ISIS 9000)
 - __: 38 (CADAS 200)
- 2. Insert indicator cuvette of the TC cuvette combination.
- 3. Insert indicator cuvette of the TIC cuvette combination.

CADAS 30S/50S, ISIS 9000:

The results for TC, TIC and TOC can be called up by pressing the key under the symbol \rightarrow .

CADAS 30S/50S/200, LASA 50/100, XION 500:

The result is displayed in TC, TIC and TOC.

Parameter	Meas. range
TOC	60 – 735 mg/L
TC	75 – 750 mg/L
TIC	15 – 150 mg/L

TOC Total organic carbon

Edition 96/03

Evaluation

- 1. Check program control number:
 - __: 38 (CADAS 200)
 - __: 34 (ISIS 6000) ⇒ Select »CUVETTE TEST« mode.
- 2. Select test number (see below).
- 3. Control number must be 2.
- Insert indicator cuvette of the TC cuvette combination and press green key.
- Insert indicator cuvette of the TIC cuvette combination and press green key.

ISIS 6000:

The results for **TC**, **TIC** and **TOC** can be called up by pressing the key under the symbol \rightarrow .

CADAS 200:

The result is displayed in TC, TIC and TOC.

Meas. range
60 – 735 mg/L
75 – 750 mg/L
15 – 150 mg/L

TOC Total organic carbon

Edition 96/03

Evaluation

- 1. Select »TEST« mode.
- 2. Select symbol (see below).
- 3. Select symbol » > «.
- 4. Check factors and measuring wavelength in memory »Mem«.
- Insert blank-value cuvette (distilled water) and press "NULL" (zero) key.
- Insert indicator cuvette of the TC cuvette combination and press "MESS" (measure) key.
- Insert blank-value cuvette (distilled water) and press "NULL" (zero) key.
- 8. Insert indicator cuvette of the **TIC** cuvette combination and press "MESS" (measure) key.

The result is printed out in TC, TIC and TOC.

If more than one sample is to be measured start the next evaluation at point 5.

Parameter	Symbol	Meas. range
TOC	\$ 381	60 – 735 mg/L
TC		75 – 750 mg/L
TIC		15 – 150 mg/L

CADAS 100 (≥ LPG 210)

LCK 381

TOC Total organic carbon

Edition 96/03

Evaluation

- 1. Select »TEST« mode.
- 2. Select symbol (see below).
- 3. Control number must be 8.
- 4. Insert blank-value cuvette (distilled water) and press "NULL" (zero) key.
- 5. Insert indicator cuvette of the **TC** cuvette combination and press "MESS" (measure) key.
- Insert indicator cuvette of the TIC cuvette combination and press "MESS" (measure) key.

The result is printed out in TC, TIC and TOC.

If more than one sample is to be measured start the next evaluation at point 5.

Parameter	Symbol	Meas. range
TOC	381	60 – 735 mg/L
TC		75 – 750 mg/L
TIC		15 – 150 mg/L

DR 1900

LCK 381

TOC Total organic carbon

Edition 96/03

Evaluation

- 1. Select »Barcode Programs«.
- 2. Select test number (see below).
- 3. Control number must be 2.
- Insert indicator cuvette of the TC cuvette combination and press »Read 1«.
- Insert indicator cuvette of the TIC cuvette combination and press »Read 2«.

The result is displayed in TC, TIC and TOC.

Parameter	Test-No.	Meas. range
TOC	381	60 – 735 mg/L
TC		75 – 750 mg/L
TIC		15 – 150 mg/L

