
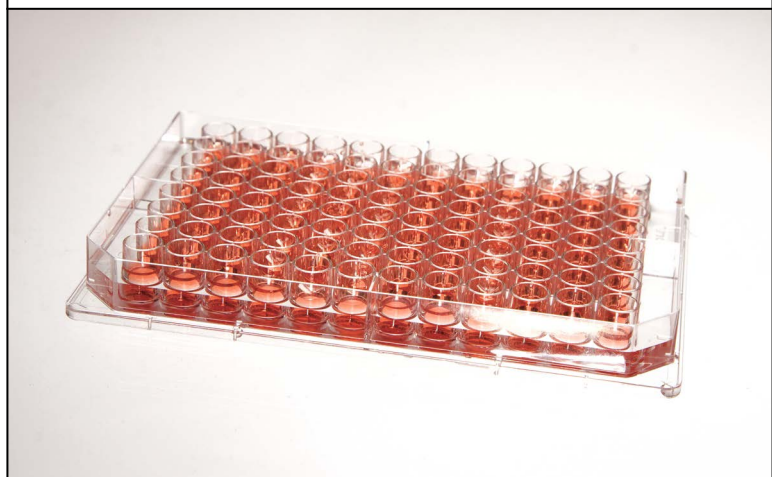
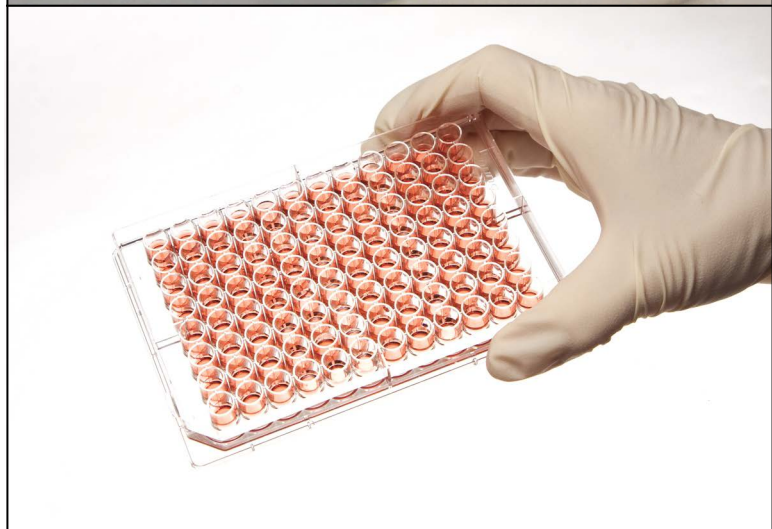


# Nunc™ Microplate Edge 2.0

up to  
**65%  
off**

When it comes to growing healthy, viable cells in extended microplate cultures, evaporation is a pressing issue that researchers are commonly faced with. In order to mitigate the damages caused by evaporation, notably among the highly vulnerable perimeter wells, researchers typically leave the outer wells in a standard 96-well plate empty ultimately sacrificing 37.5% of each plate. With the Thermo Scientific™ Nunc™ Edge 2.0 Plate, this work-around strategy is no longer necessary.

- Edge plate reduces edge effect for full plate results
- Reclaim full-use of your 96-well plate
- Maintain cell viability and consistency across all 96 wells of your plate.
- Particularly suited to long term assays where evaporation is an issue
- Reduces overall plate evaporation to <2%, yielding consistent, viable cells across the entire 96 well plate
- Using all 96 wells reduces waste
- Simple to use, just add 1.7mL of liquid to each quadrant in the outer reservoir

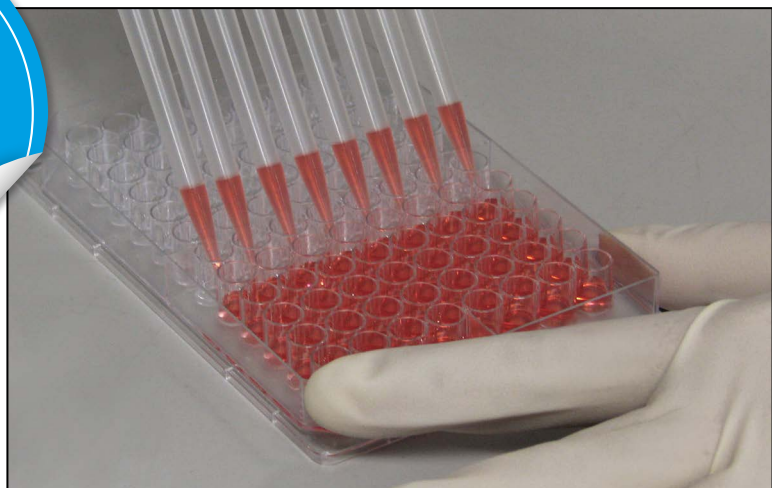

Offer valid until 30<sup>th</sup> September 2017

# Nunc™ Microplate Edge 2.0

Code	Alt Ref	Description	Pack	List Price £	Offer Price, £
WEB1700	167425	96 well, flat bottom, clear, Nunclon™ Delta, sterile, with lid, individually bagged	50	132.70	46.91
WEB1738	167542	96 well, flat bottom, clear, Nunclon™ Delta, sterile, with lid, bulk pack, 10 per bag	160	329.29	135.70
WEB1740	167554	96 well, flat bottom, clear, Nunclon™ Delta, sterile without lid, bulk pack, 10 per bag	160	306.09	126.14
WEB1742	167574	96 well, flat bottom, clear, Nunclon™ Delta, sterile, without lid, individually bagged	50	113.70	40.17
WEB1744	267427	96 well, flat bottom, clear, non-treated, sterile, with lid, individually bagged	50	105.80	37.39
WEB1746	267544	96 well, flat bottom, clear, non-treated, sterile, with lid, bulk pack, 10 per bag	160	272.94	112.48
WEB1748	267556	96 well, flat bottom, clear, non-treated, sterile, without lid, bulk pack, 10 per bag	160	254.14	104.73
WEB1750	267566	96 well, flat bottom, clear, non-treated, non sterile, without lid, bulk pack, 10 per bag	160	249.73	102.91
WEB1752	267576	96 well, flat bottom, clear, non-treated, sterile, without lid, individually bagged	50	92.50	32.58
WEB1754	267578	96 well, flat bottom, clear, non-treated, non sterile, with lid, bulk pack, 10 per bag	160	245.31	101.00

## FREE samples on request

Please contact your local Territory Sales Manager

Offer valid until 30<sup>th</sup> September 2017