


## SLS Flowgen Mini Rapide Horizontal Electrophoresis System



An ultra-compact self-contained system for routine molecular biology procedures and quick checks of samples. Buffer and gel volumes kept to a minimum to maximise current and separation speed. Parallel electrode arrangement allows ultra efficient current transfer, enabling resolution to be completed within 30 minutes. The UV transparent base allows direct viewing on a UV transilluminator with no need for time consuming transfer and potential gel damage.

- Convenient self contained solution for routine agarose gel electrophoresis
- Ultra compact
- Low buffer volumes
- Dual comb slots allow the loading of up to 40 samples per gel
- Supplied with 2 x 8 sample combs (1.5mm) and casting dams

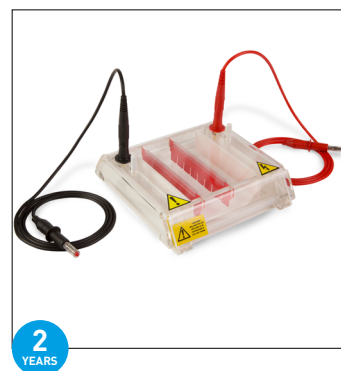

Code	Gel Dims, w x l, mm	Dims, w x d x h, mm	Max. Capacity	Base Buffer Volume, mL	Price
<b>FMMS10</b>	100 x 80	150 x 150 x 40	40 samples	50	£332.00

### Combs – 1mm Thick

Colour coded: white.

Code	Description	Price
<b>FMMS11</b>	1 sample, 1 dye	£29.00
<b>FMMS41</b>	4 sample	£29.00
<b>FMMS81</b>	8 sample	£29.00
<b>FMMS121</b>	12 sample	£29.00
<b>FMMS161</b>	16 sample	£29.00
<b>FMMS201</b>	20 sample	£25.20

### Combs – 1.5mm Thick

Colour coded: red.

Code	Description	Price
<b>FMMS115</b>	1 sample, 1 dye	£27.00
<b>FMMS415</b>	4 sample	£27.00
<b>FMMS815</b>	8 sample	£29.00
<b>FMMS1215</b>	12 sample	£29.00
<b>FMMS1615</b>	16 sample	£29.00
<b>FMMS2015</b>	20 sample	£23.10

### Accessories

Code	Description	Pack	Price
<b>FMMSDAM</b>	Casting dams	2	£36.00
<b>RPW02</b>	Replacement platinum wire, 0.2mm - 50cm	Each	£108.00