


NEW

SLS Select EasySeal qPCR Crystal Adhesive Microplate Film



- Optically clear adhesive film, pressure activated adhesive, peelable
- Suitable for qPCR and other imaging techniques including crystallisation
- When the seal is in position, pressure can be applied to burst the adhesive capsules, releasing a strong adhesive only where the seal touches the raised well rims of the plate – the rest of the seal area above the wells remains optically clear

Code	Dims, w x l, mm	Pack	Price
SLS9205	80 x 140	100	£99.00

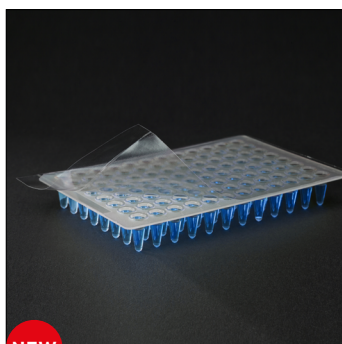


NEW

SLS Select EasySeal Micro Adhesive Microplate Film



- Low strength adhesive film, transparent and peelable, suitable for short term storage
- A low cost sealing option, suitable for temporary storage and as a cover for applications such as centrifugation

Code	Dims, w x l, mm	Pack	Price
SLS9220	80 x 135	100	£42.90

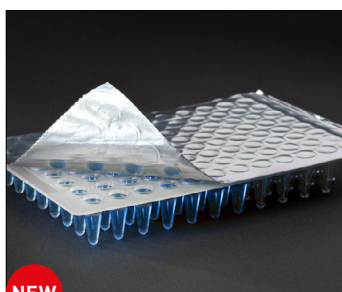


NEW

SLS Select EasySeal Foil PCR Microplate Film



- Pierceable adhesive aluminium foil, strong adhesive and peelable
- Suitable for PCR and other high temperature applications
- When pierced, the foil tears in an irregular manner which prevents the formation of a vacuum
- Perforated end tabs for easy application and removal by peeling

Code	Dims, w x l, mm	Pack	Price
SLS9238	80 x 130	100	£86.00

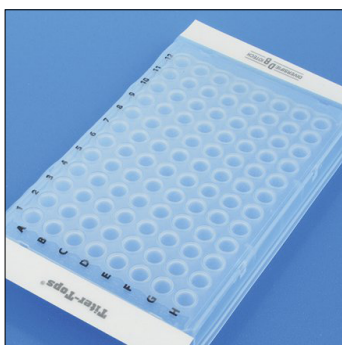

Titer-Tops® Microplate Sealing Film

- Durable and thick (4.4mil) yet easily pierced
- Prevents well-to-well contamination, sample spillage and evaporation
- Non cytotoxic, sterile, solvent resistant and waterproof polythene
- Resistant to 100% humidity and temperatures from -40 to 100°C
- Chemically resistant and thermally stable
- Easily marked with any writing instrument
- Two scored end taps for easy positioning
- UV transparent down to 250nm
- Sheet dims, w x l: 152 x 83mm

Code	Alt Ref	Pack	Price
MIC5000	T-TOPS-50	50	£36.60
MIC5002	T-TOPS-100	100	£62.30