



SLS Lab Basics DigiRoll Tube Roller



Digital tube rollers with microprocessor control for optimum mixing of viscous samples and accurate speed under variable load. The brushless DC motor for maintenance free long life and quiet operation. Stackable design with robust construction.

- Rocking and rolling motion
- Digital display with soft touch keypad
- IP 21 protection class
- Timer from 1 min to 99 min and infinite mode
- User can set pulse time from 30 to 99 seconds
- Adjustable rollers spacing to hold different diameter tubes
- User can save and run up to 9 user defined programs

Code	No. of Rollers	Capacity	Speed Range, RPM	Tilt Angle, °	Dims, w x d x h, mm	Price
SLS4234	6	15 x 3/4mL tubes	10 - 80	±3	502 x 301 x 93	£473.00




SLS Lab Basics 10 Roller Tube Roller



Designed for gentle, thorough mixing. As the roller mechanism rotates through each complete revolution, the right side of each roller is gently elevated about 13mm and returns to the horizontal starting position. This additional, gentle up/down motion provides a "wave" effect and guarantees that the entire sample is subjected to complete, thorough mixing, end to end. This built in tilt feature provides a level of homogeneous mixing uniformity that cannot be matched by any stationary roller apparatus.

- Gentle rolling motion
- Smooth ±5° tilt action
- For tubes/bottles from 8 to 200mm diameter

Code	No. of Rollers	Speed, RPM	Tilt Angle, °	Dims, w x d x h, mm	Price
SLS7048	10	38	±5	450 x 180 x 90	£877.00

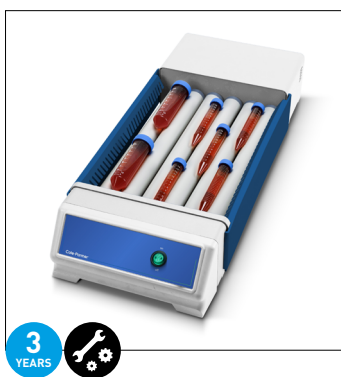

Cole-Parmer® SRT Series Rollers



Solidly constructed roller shakers, recommended for mixing samples, viscous substances and liquid solid suspensions where minimum aeration is required. Shakers accommodate most tube and bottle sizes.

- Gentle rocking and rolling action for complete mixing
- Safe for use in cold rooms down to 4°C or in incubators to 60°C, up to 80% humidity
- Continuous and quiet operation
- Choice of analog fixed speed or digital variable speed models

Code	Alt Ref	Model	No. of Rollers	Control	Speed Range, RPM	Dims, w x d x h, mm	Price
MIX1982	04653-02	SRT6	6	Analog	33 (fixed)	595 x 225 x 110	£716.00
MIX1986	04653-06	SRT6D	6	Digital	5 - 60	595 x 255 x 110	£944.00
MIX1984	04653-04	SRT9	9	Analog	33 (fixed)	595 x 375 x 110	£870.00
MIX1988	04653-08	SRT9D	9	Digital	5 - 60	595 x 375 x 110	£1095.00