



Warranty excludes rechargeable battery (1 year warranty) and replacement filter


SLS Lab Basics Pipette Fillers



Lightweight, made from UV resistant material, features adjustable valves. Filling and dispensing speed is easily set using the same hand by pressing the large mode button. The bright LEDs will indicate which of the three modes has been selected from high, low or gravity mode. User labelling area provides customisation, instant recognition and/or task designation. Additional features include hand neutral design, concave padded triggers, a 1mL pipette adapter to increase stability of small pipettes, a long life Li-Ion battery that allows 20 hours of continuous usage, and low battery indicator.

- For 1 - 100mL glass or plastic pipettes
- Uses standard 0.45µm hydrophobic filters
- Supplied with charging stand, magnetic/fixable wall mount, nozzle adapter, rechargeable battery, 0.45µm filter, 5 x coloured nose cones and low voltage double insulated power adapter

Code	Description	Price
SLS7010	Pipette filler	£200.00
PIP6042	Sterile filter, 0.45µm	£11.45

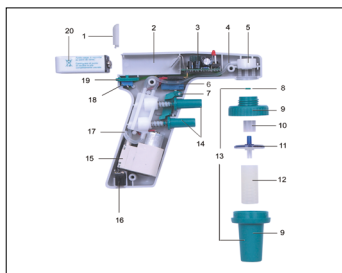

INTEGRA PIPETBOY acu 2 Pipette Fillers



For use with serological pipettes in the volume range 1 to 100mL: pipetting is controlled by two operating buttons and an adjustment wheel for setting the pipetting speed. Precise drop-by-drop dispensing is easily possible by applying a slight pressure on the operating button. The PIPETBOY acu 2 is powered by high-performance lithium polymer battery - power for more than 5500 cycles.

- Autoclavable sterile module, easy to replace filter
- LED provides information about charging status
- Continuous preset of speed via thumb wheel
- Sensitive control by applying finger pressure
- Lightweight and well-balanced
- Maximum pipetting speed: 13.5mL/s (with a 50mL pipette)
- Kit consists of PIPETBOY acu 2, rechargeable 9V lithium polymer battery, power supply, sterile 0.45µm filter, wall mount and free Swiss pocket tool

Code	Description	Price
I155000	Classic kit	£416.00
I155015	Green kit	£416.00
I155016	Red kit	£416.00
I155017	Blue kit	£416.00
I155018	Transparent kit	£416.00
I155019	Purple kit	£416.00
I155024	Pink sunset kit	£416.00

Accessories/Spare Parts

Code	Description	Price
I153010	Separate rapid charger for 2 batteries (incl. one battery)	£81.00
I155160	Battery compartment lid for PIPETBOY green	£9.00
I155161	Battery compartment lid for PIPETBOY red	£9.00
I155162	Battery compartment lid for PIPETBOY blue	£9.00
I155163	Battery compartment lid for PIPETBOY transparent	£9.00
I155203	Battery compartment lid for PIPETBOY grey	£9.00
I155164	Battery compartment lid for PIPETBOY purple	£9.00
I155150	Handle, without components for PIPETBOY green	£36.00
I155151	Handle, without components for PIPETBOY red	£36.00
I155152	Handle, without components for PIPETBOY blue	£36.00
I155153	Handle, without components for PIPETBOY transparent	£36.00
I155215	Handle, without components for PIPETBOY grey	£36.00
I155154	Handle, without components for PIPETBOY purple	£36.00