


## Syphon Pump

Syphon pump with a polyethylene squeeze bulb on a 500mm vertical tube, and a 540mm flexible side arm for dispensing from bulk. Syphons water at approx 2.7L per minute.

Code	Material	Length, mm	Price
PUM1004	Polyethylene	500	£14.12

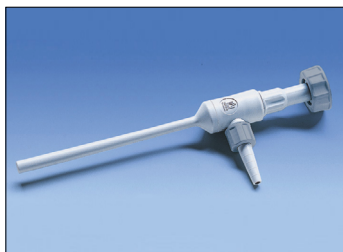

## BRAND Water Jet Pump



For generating vacuum and extracting liquids and vapours. Very low water consumption, constant ultimate pressure and high suction capacity. High chemical resistance, fluid path consists of PP, FKM and PTFE, integrated non-return valve increases safety. Detachable vacuum connection, R  $\frac{3}{4}$ " connector fitted. Operating temperature up to a maximum of 80°C.

Code	Alt Ref	Material	Length, mm	Price
FIL1110	159600	Polypropylene	210	£51.00



## Kartell Water Jet Pump



This pump will work on pressures up to 10kg/cm<sup>2</sup>, and a built in non return valve eliminates back flow. Useful for sucking to waste any split chemicals as well as for filter work. Easily dismantled for cleaning. Autoclavable.

Code	Alt Ref	External H <sub>2</sub> O Tube Ø, mm	External Vacuum Tube Ø, mm	Price
FIL1106	1395	11.85	9	£32.00



## Metal Water Pumps

Constructed in nickel plated brass and stainless steel. Designed on conventional water-jet lines to give vacuum readings of up to 625mm. Fitted with riffled nozzles for vacuum and water exit nozzles. Connection to water tap can be plain for rubber tubing connection or with threaded wing nut. Neoprene non-return valve fitted.

Code	Description	Price
PUM1200	Plain connection	£57.00
PUM1202	$\frac{1}{2}$ " BSP connection	£74.00



## SLS Lab Basics Mini Vacuum Pumps



Single-head, dry-running devices used in a wide range of laboratory applications. They transfer, compress and pump down without contamination.

- 100% oil-free transfer
- PTFE coated version for slightly aggressive or corrosive gases and vapours
- Maintenance-free
- High level of gas tightness
- IP20 motor protection
- Silencer fitted to outlet port
- 6mm hosetail on inlet

Code	Diaphragm Material	Delivery, L/min	Ultimate Vacuum, mbar	Price
SLS2600	EPDM rubber	6	100	£325.00
SLS2602	PTFE coated	5.5	160	£383.00