



Restriction Enzyme Hpy99 I




Cat.# FG-Hpy99I	Size 100 units	Conc. 2 units/μl
---------------------------	--------------------------	----------------------------

Store at -20°C

Supplied with: 10X FastGene® Buffer IV (FG-REB4)
10X FastGene® FastCut Buffer (FG-REBHF)
6X DNA Loading Buffer
Sterile water

Recognition site

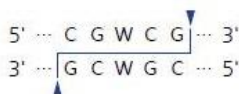

For Research Use Only. Not for use in diagnostic procedures.



Dilution buffer:

FastGene® Diluent A

Heat Inactivation

Hpy99 I can be inactivated at 65°C for 20 min.

Methylation sensitivity

dam methylation: Not sensitive

dcm methylation: Not sensitive

CpG methylation: sensitive

Prolonged incubation

A minimum amount of enzyme required to digest 1 μg substrate DNA for 16 hr; 1 U.

Relative activity in FastGene® Buffers

FastGene® Buffer I:	100%
FastGene® Buffer II:	25%
FastGene® Buffer III:	10%
FastGene® Buffer IV:	100%
FastGene® FastCut Buffer:	100%

Note

Cleavage of mammalian genomic DNA is blocked by CpG methylation.

Source: *Helicobacter pylori* J99

Reaction conditions

1X FastGene® Buffer IV 37°C

1X FastGene® FastCut Buffer, 37°C

FastGene® FastCut Buffer

FastGene® restriction enzyme can cut substrate DNA in 5-15 with FastGene® FastCut Buffer.

1X FastGene® Buffer IV

20 mM Tris-acetate (pH 7.9 at 25°C)

50 mM potassium acetate

10 mM magnesium acetate

100 μg/ml BSA

Unit definition

One unit is defined as the amount of enzyme required for complete digestion of 1 μg bacteriophage λ at 37°C for 1 hr in 50 μl reaction mixtures.

Quality control

- Unit definition assay
- Overdigestion assay
- Endonuclease assay
- Extreme pure assay

Standard reaction condition

- Normal protocol

Component	Final Conc.	Volume
Substrate DNA	1 μg	X μl
10X FastGene® Buffer IV	1 X	5 μl
Hpy99 I	2 unit	1 μl
Sterile water		up to 50 μl

→ Incubate at 37°C for 1 hr

- Fast protocol

Component	Final Conc.	Volume
Substrate DNA	1 μg	X μl
10X FastGene® FastCut Buffer	1 X	5 μl
Hpy99 I	2 unit	1 μl
Sterile water		up to 50 μl

→ Incubate at 37°C for 15 min

※ We recommend 5-10 units of enzyme per μg DNA and 10-20 units for genomic DNA in a 1 h digest.