



Moticam 1080 HDMI

FULL HD CAMERAS

The Full HD affordable multi-tasking microscopy camera. The Moticam 1080 is dedicated to documentation of microscopic results with a clear focus on maximum colour fidelity and fast live image.

It does not require a computer for standard operation, as it comes with an on-board software that will allow you to use it by simply connecting it to a screen through its HDMI port, and with a mouse connected to the USB port of the camera.

Moticam

www.moticamseries.com | www.moticeurope.com

FULL HD STANDALONE CAMERA





ON-BOARD SOFTWARE & MOTIC IMAGES PLUS 3.0 SOFTWARE

Moticam¹⁰⁸⁰ TECHNICAL SPECIFICATIONS


CMOS
1/2.8"
8.0 MP
Still Image 3840 x 2160 (8.0MP) / Video Full HD 1920 x 1080
1920 x 1080 (Full HD) @ 30 fps*
1920 x 1080 (Full HD) @ 60 fps*
2.8µm x 2.8µm
HDMI (1080p) and USB2.0
SD Card (max 32Gb)
Microsoft Windows XP SP3/Vista/7/8/10 and MAC OSX & Linux
2GHz dualcore – RAM memory 2GB – Video memory min. 512 MB
C-Mount
Standard tools Still image capture, Live image record, Zoom in & out, Magnifier, ROI, Mirror, Rotation, Split, Freeze, Cross line, Gallery.
Image adjustments (AE/AWB) Auto/Manual Exposure mode, Gain, Target Brightness setting, Light Frequency, Auto/Manual White Balance mode.
Advanced settings Saturation, Contrast, Gamma, Sharpness, Denoise, User settings, System's Time setting, Language (English, Chinese and Japanese), Format SD card.
CS Ring Adaptor, HDMI cable, Focusable Lens, 30mm and 38mm
Eyepiece Adapters, USB 2.0 Cable, Calibration Slide, Macro Tube,
Motic Images Plus 3.0 for PC/OSX/Linux


under

THE [AllinOneBOX] CONCEPT Everything that you may need to work with your Moticarm is included in one box

On-board software

The Moticam 1080 has its own built-in software, making it possible to view live images, make image adjustments or edit them. You can also capture images or record videos directly from a screen with the help of a mouse, storing them in an SD card.

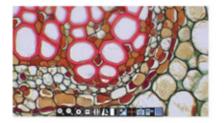

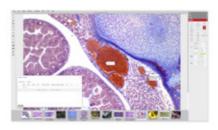

Motic Images Plus 3.0

By connecting your Moticam to your computer you will be able to work with our well-known software, that comes included in the package. View, capture, edit, measure, make reports... all its standard features and the new ones packed in a new user-friendly interface.