Axygen® Power Supplies

PS300GP and PS250HA Ideal for most of your electrophoresis needs



The Axygen Power Supplies are full-featured power sources ideal for molecular biology workflows including DNA, RNA, and protein gel electrophoresis applications, as well as blotting. The compact user interface has 4 sets of color-coded terminals allowing for use with up to 4 gel boxes simultaneously. Depending on the model, operating modes include constant voltage, current, or power all with automatic crossover feature. The Power Supplies feature an easy to set, LCD backlit, digital display to perform according to the needs of your operational parameters. The PS250HA has the ability to set programs with 10 steps per program. Both units have timers, are stackable, and are engineered with safety in mind to meet the requirements that todays protocols demand.

PS250HA

- High amperage great for Western blotting
- ▶ Programmable

PS300GP

- Ideal for small format protein gels or DNA horizontal gels
- Easy to run constant mode for consistent Voltage or Amperage protocols



Specifications

	Axygen PS250HA Power Supply	Axygen PS300GP Power Supply
Input Power (switchable)	115 VAC, 50/60 Hz	115 VAC, 50/60 Hz
	230 VAC, 50/60 Hz	230 VAC, 50/60 Hz
Output power	300W	90W
Output voltage range	5 to 250V	2 to 300V
Output current range	10 to 3000 mA	4 to 500 mA
Timer	Approx. 99.99 hr./min.	Approx. 99.99 hr./min.
Terminal pairs	4 positive voltage and 4 negative voltage	4 positive voltage and 4 negative voltage
Operating Modes	Constant Voltage: 1V step Constant Current: 10 mA step Constant Power: 1W step	Constant Voltage: 1V increment settings Constant Current: 1 mA increment settings
Crossover	Automatic	Automatic
Pause function	Yes	Yes
Safety features	No Load Detection Load Change Detection Overload Detection Ground Leak Detection Auto Restart	No Load Detection Load Change Detection Overload Detection Ground Leak Detection Auto Restart
Programmable	Yes. Can store up to ten 10-step methods	No
Housing size (W x D x H)	7.9 x 12.4 x 3.4 in. (20 x 31.4 x 8.5 cm)	7.9 x 12.4 x 3.4 in. (20 x 31.4 x 8.5 cm)
Weight	5.4 lbs (2.45 kg)	5.4 lbs (2.45 kg)
Certifications	CE; cETLus	CE; cETLus
Warranty	3 years	3 years



Ordering Information

Axygen® Power Supplies

Cat. No.	Description	
PS250HA	Axygen 250V High Amperage Power Supply, 115V	
PS250HA-230V	Axygen 250V High Amperage Power Supply, 230V	
PS300GP	Axygen 300V General Purpose Power Supply, 115V	
PS300GP-230V	Axygen 300V General Purpose Power Supply, 230V	


*NOTE: The 115V and 230V in the product description refer to the plug type necessary to operate in your region.

Related Products

Cat. No.	Description	
HGB-7	Axygen Horizontal Gel Box (7 cm)	
HGB-10	Axygen Horizontal Gel Box (10 cm)	
HGB-15	Axygen Horizontal Gel Box (15 cm)	
HGB-20	Axygen Horizontal Gel Box (20 cm)	
GD-1000	Axygen Gel Documentation System	
GDBL-1000	Axygen Gel Documentation System with Epi blue light	





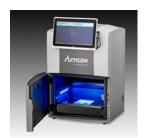

For more specific information on claims, visit www.corning.com/certificates.

Warranty/Disclaimer: Unless otherwise specified, all products are for research use only. Not intended for use in diagnostic or therapeutic procedures. Corning Life Sciences makes no claims regarding the performance of these products for clinical or diagnostic applications.

For additional product or technical information, visit www.corning.com/lifesciences or call 800.492.1110. Outside the United States, call +1.978.442.2200 or contact your local Corning sales office.

CORNING

Corning Incorporated Life Sciences

836 North St.
Building 300, Suite 3401
Tewksbury, MA 01876
t 800.492.1110
t 978.442.2200
f 978.442.2476
www.corning.com/lifesciences

ASIA/PACIFIC Australia/New Zealand t61 427286832 Chinese Mainland t86 21 3338 4338 f 86 21 3338 4300

India t 91 124 4604000 f 91 124 4604099 Japan t 81 3-3586 1996 f 81 3-3586 1291 Korea t 82 2-796-9500

f 82 2-796-9300 **Singapore** t 65 6572-9740 f 65 6735-2913

Taiwan t 886 2-2716-0338 f 886 2-2516-7500 EUROPE CSEurope@corning.com

t 0800 916 882 f 0800 918 636 **Germany** t 0800 101 1153

France

f 0800 101 2427

The Netherlands
t 020 655 79 28
f 020 659 76 73

United Kingdom

United Kingdom t 0800 376 8660 f 0800 279 1117 All Other European Countries t +31 (0) 206 59 60 51

t+31 (0) 206 59 60 51 f+31 (0) 206 59 76 73

LATIN AMERICA grupoLA@corning.com Brazil t 55 (11) 3089-7400 Mexico t (52-81) 8158-8400