



DATA SHEET

Product Code SLS1242 – SLS Select 1.2mL round well 'U' bottom plate no sealing rim


About SLS1242

The 1.2mL round well 'U' bottom plate has a total volume of 1.24mL. The top of the plate is flat without heat sealing rims. SLS Select also provides a compatible, reusable silicone cap mat that can be punctured (SLS1288).

The plate is manufactured in ISLS Select's class 8 cleanrooms using a medical-grade polypropylene. This polypropylene helps achieve very low concentrations of extractable elements and has a temperature capability range of -80°C to +121°C.

This plate has been developed precisely to conform to the footprint dimensions defined by SLAS/ANSI standards. This guarantees its compatibility with both manual and automated sample handling systems, readers, and microplate washers.

SLS Select also offer a version of this plate with rims. Product code SLS1240.

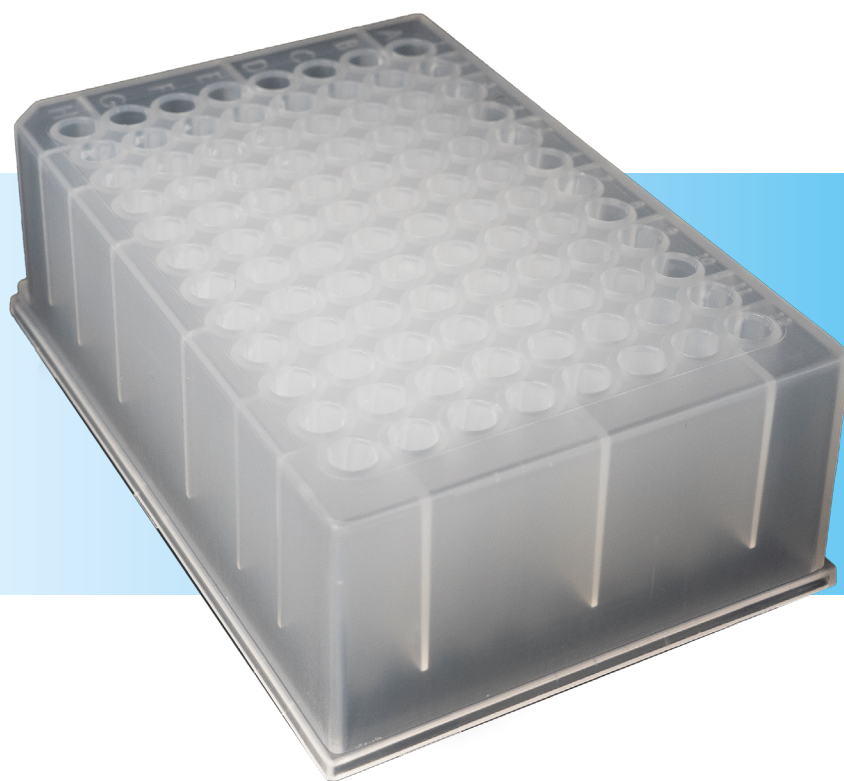

United Kingdom:
slsinfo@scientific-labs.com

International:
slsexport@scientific-labs.com

Northern Ireland:
NI.sales@scientific-labs.com

Republic of Ireland:
ROI.sales@scientific-labs.com

 **SCIENTIFIC
LABORATORY
SUPPLIES**



DATA SHEET

Product Code SLS1242 – SLS Select 1.2mL round well 'U' bottom plate no sealing rim

KEY FACTS:

- Manufactured in ISO Class 8 Cleanroom
- Alphanumeric well coding
- High temperature and chemical resistance (-80°C to +121°C)
- Autoclavable
- RNase, DNase Free and Non-Pyrogenic
- Stackable for use in plate & automation hotels
- Stack Height for 2 plates 81.50mm
- Available sterile or non-sterile
- Flat surface for sealing with adhesive and heat seals




ORDERING INFORMATION:

- Case Dimensions: 47 x 47 x 18cm
- Case Weight: 7kg
- Case qty: 50 Plates
- Packsize: 10 Bags x 5 Plates
- Standard pallet quantity: 34 cases

SLS code	Qty per pack	Wells	Well bottom	Max volume mL	L	W	H	Well depth	Well diameter	Stack height of 2 plates	Silicone cap matt
SLS1242	50	96	U	1.24	127.75	85.48	42	38.5	7	81.5	SLS1288

PRODUCT DRAWING: SLS1242

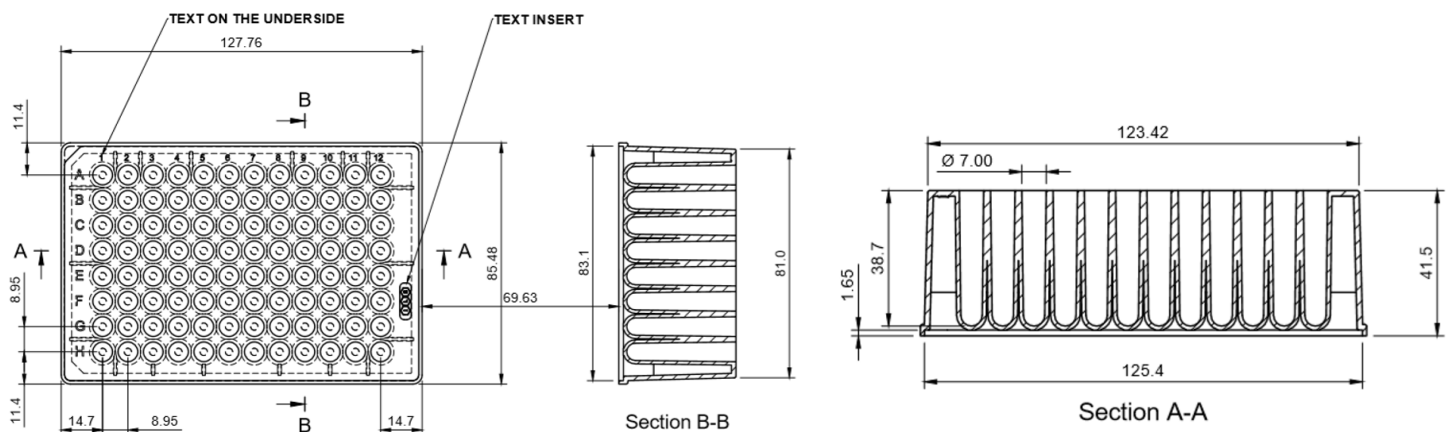

United Kingdom:
slsinfo@scientific-labs.com

International:
slsexport@scientific-labs.com

Northern Ireland:
NI.sales@scientific-labs.com

Republic of Ireland:
ROI.sales@scientific-labs.com

